



LEONE CLASS II CORRECTOR

THE NEW NON-COMPLIANCE DEVICE FOR TREATING CLASS II CASES

The reduction in size and the inner spring provide an optimal patient comfort




Courtesy of Dr. Ivano Matoni

There are 3 sizes available with the possibility of either mesial or distal attachment to the upper molar tube, for a total of 6 positions in the mouth.

When utilized in the mesial position, it can be anchored to a direct bonding single tube, without the need for a band or a round tube.


The mandibular advancement is physiologic and it is determined by a continuous light force delivered by the **MEMORIA®** spring of 200 g.

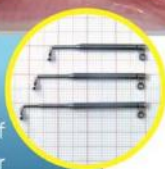

ASSEMBLY PROCEDURE WITH SUPPORT MESIAL TO VESTIBULAR TUBE



ASSEMBLY PROCEDURE WITH ANCHORAGE DISTAL TO VESTIBULAR TUBE

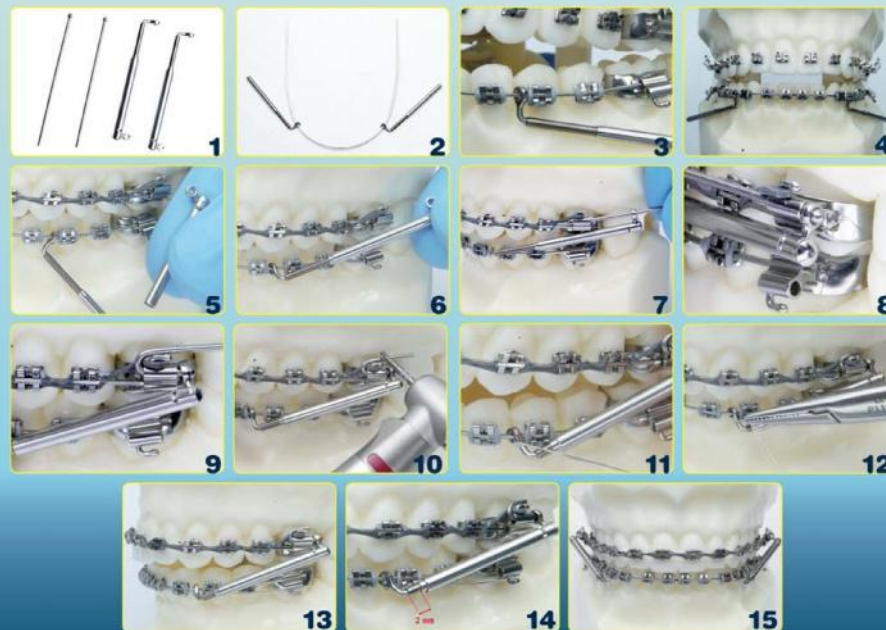

The packages include all the needful parts for the application of a bilateral corrector					
	Type	Length at maximum extension (passive configuration)	Length at maximum spring compression (without 2 mm of further possible compression)	Maximum spring extension	Maximum opening allowed to prevent disassembly
M2300-00	Extra-short	31mm	25mm	7,5mm	48mm *
M2300-01	Short	35mm	28mm	10mm	56mm *
M2300-02	Long	40mm	31mm	10mm	62mm *

*Distance between the 2 anchorage rings



LEONE S.p.a. Orthodontics and Implantology
Via P. a Quaracchi, 50 - 50019 Sesto Fiorentino
FIRENZE (Italy) - info@leone.it www.leone.it

Export sales dept:

Phone: +39.055.3044620 Fax: +39.055.304405

E-mail: export@leone.it